



ANDOVER HEALTHCARE



Wound Care Products

TABLE OF CONTENTS

COMPANY PROFILE.....	3	COFLEX TLC DRY.....	8-9
HCPCS CODES.....	4	COFLEX TLC LONG STRETCH.....	10-11
VISUAL INDICATORS.....	5	COFLEX TLC ZINC & CALAMINE.....	12-13
PRODUCT SELECTION GUIDE.....	6	COFLEX AFD.....	14
COMPRESSION GUIDELINES.....	7	COFLEX UMAFD.....	15

COMPANY PROFILE

Andover Healthcare, a leading US manufacturer of cohesive bandages and compression kits, is now part of Milliken Healthcare. The combination of the two companies gives clinicians a broader suite of products along the continuum of patient care. Together, we bring unique research and development resources: Milliken with unparalleled textile and Active Fluid Management® expertise; Andover Healthcare with leading compression capabilities. With the courage to think big, we look forward to delivering healthcare innovations that will continue to positively impact the world around us and, most importantly, patient lives.



Two Layer Compression That Meets Every Need:

- ✓ Replaces ALL 2, 3, & 4 layer kits
- ✓ Indicators on every roll
- ✓ Odor/Itch Control on Layer 1 (dry)
- ✓ Easy Hand Tear Layer 2 (all)
- ✓ Long Stretch available for inactive patients
- ✓ Same application technique across entire product line
- ✓ All products deliver consistent therapeutic compression
- ✓ Patented Zinc & Calamine Two Layer Compression



HCPCS Codes

HEALTHCARE COMMON PROCEDURE CODING SYSTEM (HCPCS) CODES



Introducing...

NO LATEX




Visual Indicators

Now on ALL CoFlex TLC Products

- ✓ Easy Application
- ✓ Intended continuous levels of compression

Take the guesswork out of application and hit the target every time!

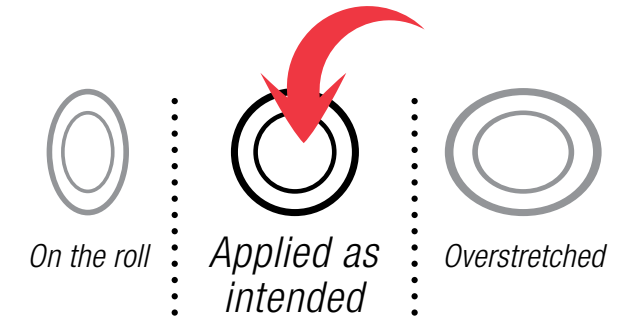

- Use visual indicators as a guide
- Stretch the bandage until the oval indicators become circles
- Visual indicators will ensure proper levels of compression based on the patient's ABI's



Product	Size	HCPCS Code
CoFlex® AFD 2"	2.5 yds.	A6209 & A6453
CoFlex® AFD 4"	2.5 yds. 7.5 yds.	A6210 & A6454
CoFlex® TLC Lite Kit CoFlex® TLC Kit	Layer 1 & Layer 2 3.4 yds. 5.1 yds.	A6441 & A6452
CoFlex® TLC XL Kit	Layer 1 & Layer 2 5.4 yds. 7 yds.	A6441 & A6452
CoFlex® TLC Zinc Lite Kit CoFlex® TLC Zinc-S Kit CoFlex® TLC Calamine Lite Kit CoFlex® TLC Calamine-S Kit	Layer 1 & Layer 2 6 yds. 7 yds.	A6456 & A6452
CoFlex® TLC LS Lite Kit CoFlex® TLC LS Kit	Layer 1 & Layer 2 4.4 yds. 8 yds.	A6441 & A6452

Centers for Medicare & Medicaid Services (CMS) has verified that the HCPCS codes listed above have been assigned to Andover Healthcare, Inc.
Updated July 2015


2 Layer Compression

PRODUCT SELECTION GUIDE

Andover Healthcare offers the widest range of two layer compression options for patients with edema associated with venous insufficiency. The systems provide a two-step level of comfort and healing unlike any other Two Layer Compression or Unna Boot on the market. All options include visual indicators for easy application.

Choose from:

- Dry or Moist (impregnated with Zinc or Calamine) comfortable foam layer
- Standard or Lite compression
- Short Stretch or Long Stretch compression

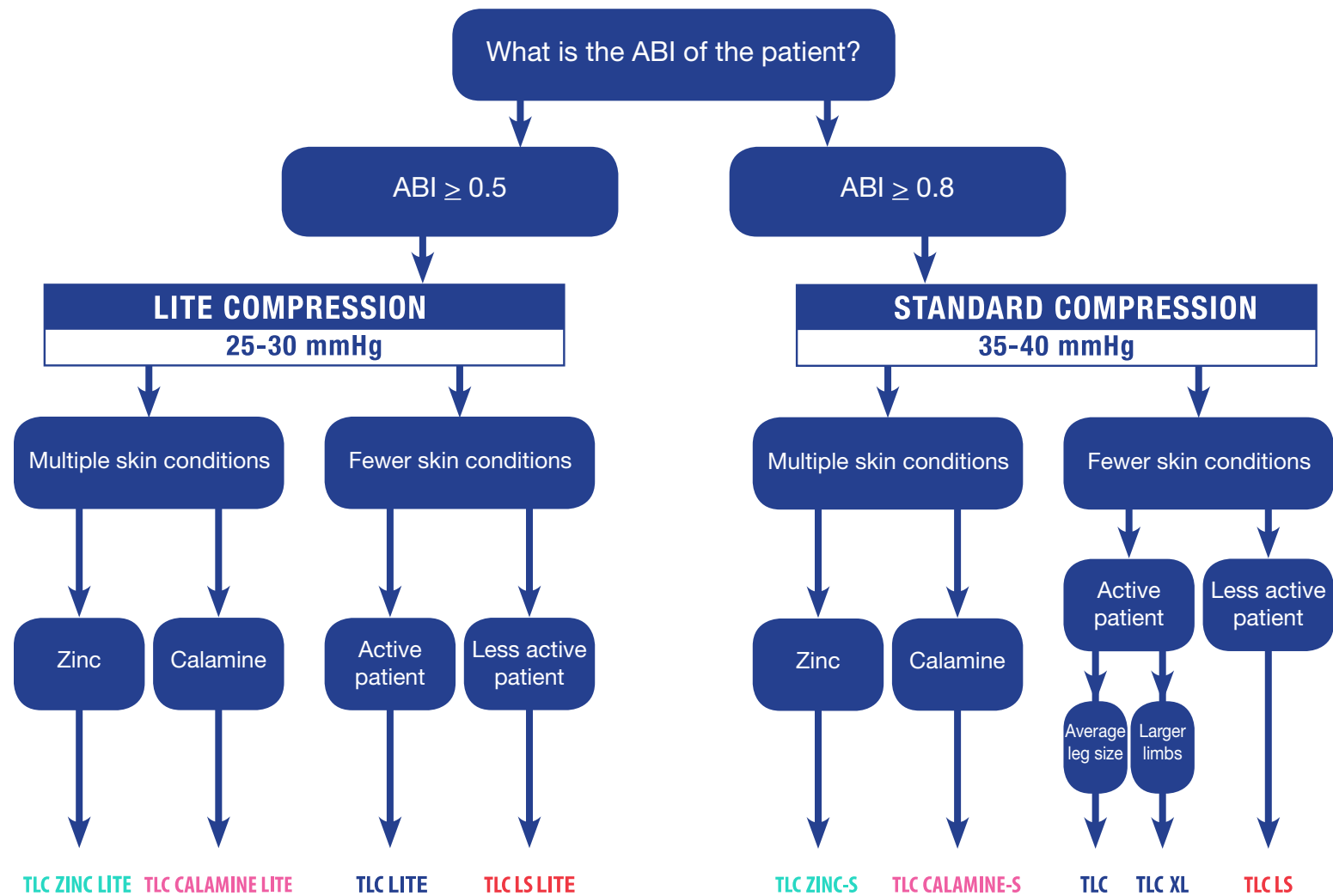

		LITE COMPRESSION 25-30 mmHg	STANDARD COMPRESSION 35-40 mmHg
SHORT STRETCH	DRY	CoFlex·TLC LITE	DRY CoFlex·TLC CoFlex·TLC XL
	MOIST	CoFlex·TLC ZINC LITE CoFlex·TLC CALAMINE LITE	MOIST CoFlex·TLC ZINC-S CoFlex·TLC CALAMINE-S
LONG STRETCH	DRY	CoFlex·TLC LS LITE	DRY CoFlex·TLC LS



CoFlex·TLC

2 LAYER COMPRESSION GUIDELINES

This product selection guide is an informational resource to explain the 9 models of CoFlex TLC and is not intended or implied to be a substitute for medical advice, diagnosis or treatment. All patients should consult their healthcare provider before making any decisions regarding compression therapy.

CoFlex® TLC

TWO LAYER COMPRESSION

DRY

NO LATEX

Visual Indicators!

Two Layer Short Stretch

CoFlex® TLC delivers a two-step short stretch performance with high working pressure and low resting pressure. Visual indicators are provided for ease of application - ovals become circles when the intended compression is achieved.

Layer 1: Soft foam with odor and itch control properties. Excellent moisture transfer rate that wicks away drainage to help keep the peri-wound skin dry.

Layer 2: Non-latex short stretch compression bandage. Sticks to itself, Easy HandTear® Technology eliminating the need for scissors.



U.S. Patent No. 5,762,523, 6,156,424
European Patent No. 0839222
European & Canadian Patent Pending

KEY FEATURES:

- Visual Indicators
- Itch Control
- Odor Control
- Easy Hand-Tear
- Nylon Stocking Included

Together the layers form a mechanical bond preventing slippage.

Use on active patients with edema.

Advantages:

- ✓ Easy to apply
- ✓ Comfortable & breathable
- ✓ Lightweight / low profile
- ✓ Sustained compression
- ✓ Non-slip support

Nylon Stocking Included

Apply nylon stocking over the completed dressing for patient comfort and ease-of-movement under clothes and on bed sheets



Available in Lite, Standard, and Standard XL.

CoFlex® TLC LITE

with Indicators

Each Box Contains 2 Rolls:
One of Each Layer + Nylon Stocking

LITE COMPRESSION
25-30 mmHg

For patients with an ABI of ≥ 0.5



Catalog #	Color	Rolls/Box	Boxes/Case	Roll Name	Description	Roll Size
7802TLC-TN	● Tan	2	8	Layer 1	Absorptive Foam Layer	4" x 3.4 yds. (Unstretched)
				Layer 2	Compression Layer	4" x 5.1 yds. (Stretched)

CoFlex® TLC

with Indicators

Each Box Contains 2 Rolls:
One of Each Layer + Nylon Stocking

STANDARD COMPRESSION
35-40 mmHg

For patients with an ABI of ≥ 0.8



Catalog #	Color	Rolls/Box	Boxes/Case	Roll Name	Description	Roll Size
7800TLC-TN	● Tan	2	8	Layer 1	Absorptive Foam Layer	4" x 3.4 yds. (Unstretched)
				Layer 2	Compression Layer	4" x 5.1 yds. (Stretched)

CoFlex® TLC XL

with Indicators

Each Box Contains 2 Rolls:
One of Each Layer + Nylon Stocking

STANDARD COMPRESSION
35-40 mmHg

For patients with an ABI of ≥ 0.8



Catalog #	Color	Rolls/Box	Boxes/Case	Roll Name	Description	Roll Size
7800TLC-XL	● Tan	2	8	Layer 1	Absorptive Foam Layer	4" x 5.4 yds. (Unstretched)
				Layer 2	Compression Layer	4" x 7 yds. (Stretched)

TLC XL has the same properties as standard CoFlex® TLC, but has an additional 2 yards on each roll for larger limbs. TLC Lite has the same first layer as standard CoFlex® TLC with a lighter weight second layer compression bandage. The lower level of compression is designed for patients with an ABI ≥ 0.5 or those who are less tolerant of compression therapy.

CoFlex® TLC LS

TWO LAYER LONG STRETCH COMPRESSION

DRY

NO LATEX

Visual Indicators!

Two Layer Long Stretch

CoFlex® TLC LS delivers a two-step long stretch performance with low working pressure and high resting pressure. Visual indicators are provided for ease of application - ovals become circles when the intended compression is achieved. Available in Standard and Lite.

Layer 1: Soft foam with odor and itch control properties. Excellent moisture transfer rate that wicks away drainage to help keep the peri-wound skin dry.

Layer 2: Non-latex long stretch compression bandage. Sticks to itself, Easy HandTear® Technology eliminating the need for scissors.



U.S. Patent No. 5,762,523, 6,156,424
European Patent No. 0839222
European & Canadian Patent Pending

Together the layers form a mechanical bond preventing slippage.

Use on less active patients with edema or where you have been using three or four layer kits.

- Advantages:**
- ✓ Easy to apply
 - ✓ Comfortable & breathable
 - ✓ Lightweight / low profile
 - ✓ Restorative compression
 - ✓ Non-slip support

Nylon Stocking Included

Apply nylon stocking over the completed dressing for patient comfort and ease-of-movement under clothes and on bed sheets



CoFlex® TLC LS LITE

with Indicators

Each Box Contains 2 Rolls:
One of Each Layer + Nylon Stocking

LITE COMPRESSION
25-30 mmHg

For patients with an ABI of ≥ 0.5

Catalog #	Color	Rolls/Box	Boxes/Case	Roll Name	Description	Roll Size
7902TLC-TN	● Tan	2	8	Layer 1	Absorptive Foam Layer	4" x 4.4 yds. (Unstretched)
				Layer 2	Compression Layer	4" x 8 yds. (Stretched)



CoFlex® TLC LS

with Indicators

Each Box Contains 2 Rolls:
One of Each Layer + Nylon Stocking

STANDARD COMPRESSION
35-40 mmHg

For patients with an ABI of ≥ 0.8

Catalog #	Color	Rolls/Box	Boxes/Case	Roll Name	Description	Roll Size
7900TLC-TN	● Tan	2	8	Layer 1	Absorptive Foam Layer	4" x 4.4 yds. (Unstretched)
				Layer 2	Compression Layer	4" x 8 yds. (Stretched)

TLC LS Lite has the same first layer as standard CoFlex® TLC LS with a lighter weight second layer compression bandage. The lower level of compression is designed for patients with an ABI ≥ 0.5 or those who are less tolerant of compression therapy.

TWO STEP EASY APPLICATION



CoFlex® TLC

TWO LAYER COMPRESSION

ZINC & CALAMINE

NO LATEX

Visual Indicators!

Two Layer Short Stretch with Zinc or Calamine

CoFlex® TLC Zinc & Calamine delivers a two-step short stretch performance with calming benefits of zinc or calamine. This unique system combines the benefits of a traditional Unna Boot with compression therapy in a much cleaner two-step process. Visual indicators are provided for ease of application - ovals become circles when the intended compression is achieved.

Layer 1: Soft foam impregnated with zinc or calamine. Soothes and calms skin with multiple wounds or other skin conditions.

Layer 2: Non-latex short stretch compression bandage. Sticks to itself, Easy HandTear® Technology eliminating the need for scissors.

Together the layers form a mechanical bond preventing slippage.



U.S. Patent No. 5,762,523, 6,156,424, 8,779,230 B2
European Patent No. 0839222
European & Canadian Patent Pending

KEY FEATURES:

- Visual Indicators
- Soothing Zinc or Calamine
- No Unna's Boot Mess
- Moisture Resistant Controlled Compression
- Easy Hand-Tear
- Nylon Stocking Included

Use on patients with edema and multiple skin conditions.

Advantages:

- ✓ Easy to apply
- ✓ Comfortable & breathable
- ✓ Absorbs 20x's its dry weight vs. traditional Unna Boot
- ✓ 50% more active ingredients than traditional Unna Boot
- ✓ Non-slip support

Nylon Stocking Included

Apply nylon stocking over the completed dressing for patient comfort and ease-of-movement under clothes and on bed sheets



Available in Lite and Standard compressions.

CoFlex® TLC ZINC with Indicators



	Catalog #	Color	Boxes/Case	Roll Name	Description	Roll Size
LITE Compression <i>For patients with an ABI of ≥ 0.5</i>	8840UBZ-TN	● Tan	8	Layer 1	Zinc Foam Layer	4" x 6 yds. (Unstretched)
				Layer 2	Compression Layer	4" x 7 yds. (Stretched)
Available in 3" = 8830UBZ-TN						
STANDARD Compression <i>For patients with an ABI of ≥ 0.8</i>	8840UBZ-SC	● Tan	8	Layer 1	Zinc Foam Layer	4" x 6 yds. (Unstretched)
				Layer 2	Compression Layer	4" x 7 yds. (Stretched)
Available in 3" = 8830UBZ-SC						

CoFlex® TLC CALAMINE with Indicators



Calamine is known for reducing itch with its soothing & cooling effect

	Catalog #	Color	Boxes/Case	Roll Name	Description	Roll Size
LITE Compression <i>For patients with an ABI of ≥ 0.5</i>	8840UBC-TN	● Tan	8	Layer 1	Calamine Foam Layer	4" x 6 yds. (Unstretched)
				Layer 2	Compression Layer	4" x 7 yds. (Stretched)
Available in 3" = 8830UBC-TN						
STANDARD Compression <i>For patients with an ABI of ≥ 0.8</i>	8840UBC-SC	● Tan	8	Layer 1	Calamine Foam Layer	4" x 6 yds. (Unstretched)
				Layer 2	Compression Layer	4" x 7 yds. (Stretched)
Available in 3" = 8830UBC-SC						

“ We wrapped a patient with the TLC Calamine - he returned today and gave us a “thumbs up.” He said that for the first time it felt great and he didn't experience any pain. His leg looks wonderful, with wound improvement and decreased edema. We're very pleased - this patient has been challenging.

Dr. Ravi K. Kamepalli, MD
Regional Infectious Disease Center
Lima, OH

“ My patients are pleased with the TLC Zinc and Calamine, especially when they see how much edema is reduced. They're easy to apply and provide good compression - they give immediate comfort.

Dr. Joseph Snow, MD, Wound Care Specialist
Concord Hospital, Wound Healing Center
Concord, NH

“ Our patients and staff love the CoFlex Unna Boot! Patients tell us CoFlex is comfortable. They don't slide or feel slimy like other products they have used. Our nurses love it for the ease and efficiency in application.

Julie Darmon, Nurse Manager
Barnes Jewish Hospital WCC
St. Louis, MO

Each Box Contains 2 Rolls - One of Each Layer + Nylon Stocking



All-in-One Foam Pad Dressing & Bandage

CoFlex® AFD's highly absorbent foam pad is attached to soft, quick-stick CoFlex® LF² for easy application. Between the LF² layer and the absorbent foam is a waterproof film layer that prevents leakage, drying out and contamination. Also available sterile.

Key Features:

- ✓ Quick & easy to apply
- ✓ High absorption and retention of fluids (equivalent to 29 gauze pads)
- ✓ Flexible, controlled compression
- ✓ Made in the USA



PRODUCT APPLICATIONS:

- Trauma Bandaging
- Post-op cranial surgery
- General wound dressing
- Cover/protect post surgery incisions
- Post bunion surgery

Catalog #	Size	Rolls/Case	Catalog #	Size	Rolls/Case	Colors
9720TN-016	2" x 2.5 yds.	16	9760TN-008	6" x 2.5 yds.	8	●
9720(S)						
9740TN-008	4" x 2.5 yds.	8	9840TN-006	4" x 7.5 yds.	6	●
9740(S)						

Choose From: Tan (TN), Sterile Tan (S)

When ordering, please indicate catalog # and color code
For example, 9720TN-016 (XX = color code)
Most products are available in bulk packaging, call for details.

All-in-One Trauma/Combat Bandage

CoFlex® UMAFD is a strong bandage with foam pad attached, plus an additional foam pad in the core for extra absorbency. It offers quick, easy application in high stress situations. It is a great alternative to existing combat bandages. Additionally, CoFlex® UMAFD has been heat tested at 130° > 90 days.

Key Features:

- ✓ Quick & easy to apply
- ✓ High absorption and retention of fluids (equivalent to 29 gauze pads)
- ✓ Flexible, controlled compression
- ✓ Able to withstand high temperatures
- ✓ Made in the USA



Catalog #	Size	Rolls/Case	Colors
9340UM-018	4" x 3 yds.	18	●
9340UM-100	4" x 3 yds.	100	●

Available in Black only

Durable trauma bandage for situations requiring additional absorption, strength for compression, and high heat/water resistance for any of the following:

- Military - Medics, Soldiers
- ER / Trauma Centers
- First Responders
- Commercial & Industrial Sites
- Emergency Medical Services
- First-Aid Kits



We were very happy with how easy the CoFlex AFD foam pad dressing was to apply. Patients we sent home with the product were able to change the dressings without difficulty and therefore more compliant. We used the CoFlex AFD foam pad dressing on many types of wounds and were pleased with the healing we saw with them all.

Amy Mahoney, Especialista em Tratamento de Feridas
Boone Hospital, Wound Clinic
Columbia, MO



MILLIKEN HEALTHCARE

9 FANARAS DRIVE, SALISBURY, MA 01952

PHONE: 978.465.0044 TOLL FREE: 800.432.6686

FAX: 978.462.0003

www.andoverhealthcare.com / healthcare.milliken.com